


Tranquil Canal Open Ear

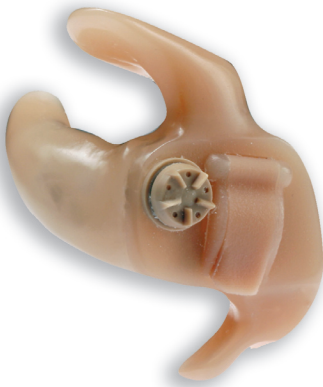

ReadyWear
Sound Generator
COE

DEVICE FEATURES

- BROAD-BAND NOISE
- #10 BATTERY
- THREE OUTPUT LEVELS AVAILABLE
- MANUAL VC
- NO OCCLUSION EFFECT
- 0 DB FLOOR LEVEL
- LOW OR HIGH FREQ. RESPONSE

CONSUMER APPROVAL

- COSMETIC APPEAL
- 1 YR. WARRANTY
- COMFORTABLE

For patients tinnitus, the COE has been the solution since 1988.

The Tranquil™ Canal Open Ear or COE, is a custom, in-the-canal custom tinnitus sound generator. The COE provides the same quality broad-band noise featured throughout the Tranquil product line while offering maximum cosmetic appeal. The free-field earmold leaves the ear canal completely open, delivering clean, stable sound while allowing wearers to hear surrounding sounds in a normal manner. Since the COE occupies only fifty percent of the ear canal, it allows for a comfortable, unoccluded fit. COE comes standard with #10 battery and 0 dB Floor volume control (no noise generator for the first 5 degrees of rotation).

When ordering the COE, specify desired output level and frequency response.

NM Circuit (normal-to-mild)

RMS Out..... 63 dB SPL
Peak Output..... 2,800 & 6000 Hz.
Battery Drain4 mA


MM Circuit (mild-to-moderate)

RMS Out..... 77 dB SPL
Peak Output..... 3,000 Hz.
Battery Drain4 mA


MS Circuit (moderate-to-severe)

RMS Out..... 94 dB SPL
Peak Output..... 3,000 Hz.
Battery Drain5 mA


NM Circuit (normal-to-mild)

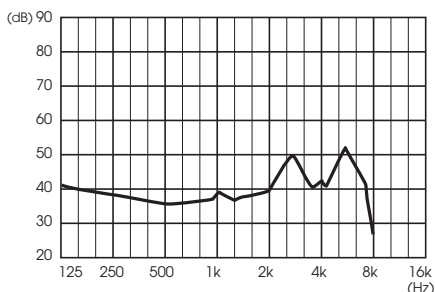

MM Circuit (mild-to-moderate)

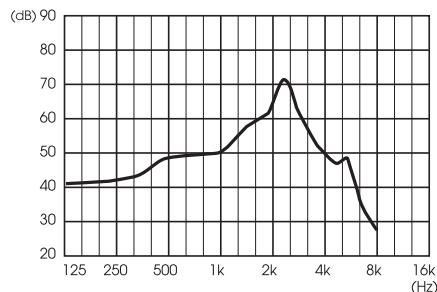

MS Circuit (moderate-to-severe)

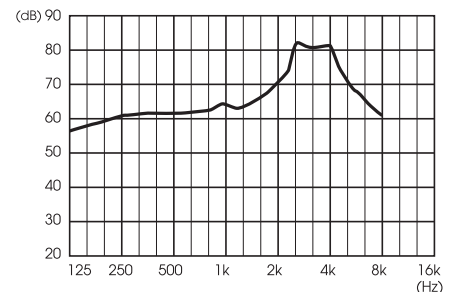

General Hearing
Instruments, Inc.

P.O. Box 23748
New Orleans, LA 70183-0748
800.824-3021 Fax (504) 733-3767